


Instructions for Rife Bioresonator


Turn the power switch to the on position.

Then Main menu appears on the display:

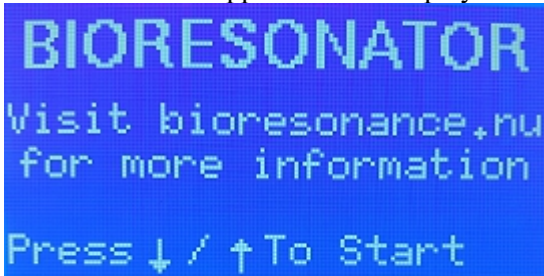

```
BIORESONATOR
Visit bioresonance.nu
for more information
Press ↓ / ↑ To Start
```

Press the up arrow to access this menu:

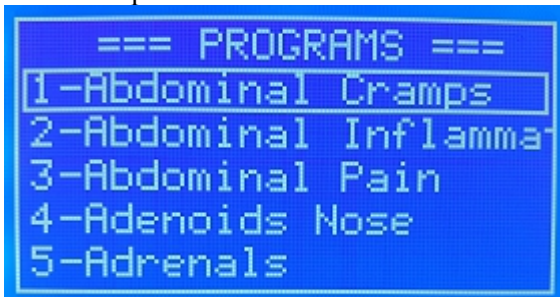


```
=== PROGRAMS ===
1-Abdominal Cramps
2-Abdominal Inflammation
3-Abdominal Pain
4-Adenoids Nose
5-Adrenals
```

In this menu it is only possible to go up or down one step at a time.
use the down and up arrow to select which program to use.


When the program is selected, (In this case program 1- Abdominal Cramps)
press start to run the selected program.

This menu will then appear:

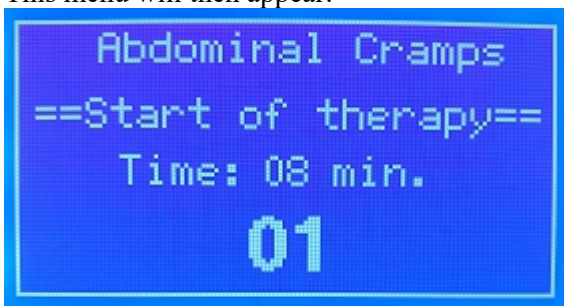

```
Abdominal Cramps
==Start of therapy==
Time: 08 min.
01
```

When the therapy program is finished, this menu appears:

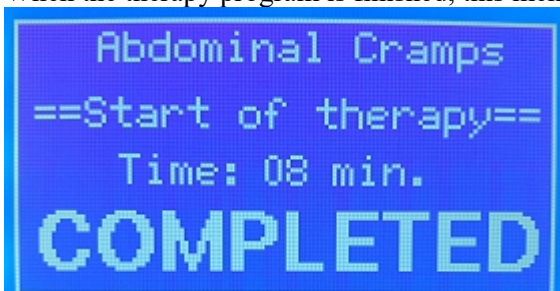

```
Abdominal Cramps
==Start of therapy==
Time: 08 min.
COMPLETED
```


Press Reset to return to main menu.

There are more than 400 pre-programmed therapy programs.
So to be able to browse faster follow the instructions below.



BIORESONATOR
Visit bioresonance.nu
for more information
Press ↓ / ↑ To Start

Press the down arrow to access this menu:

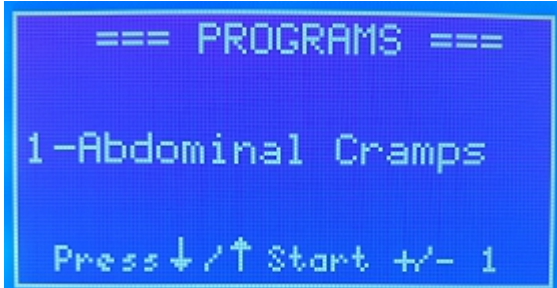

=== PROGRAMS ===
1-Abdominal Cramps
Press ↓ / ↑ Start +/- 1

Press start to choose how many steps you want to jump.

Here we have chosen 1 step at a time. Press up or down arrow to select program.



=== PROGRAMS ===
1-Abdominal Cramps
Press ↓ / ↑ Start +/- 10

Here we have chosen 10 steps at a time. Press up or down arrow to select program.



=== PROGRAMS ===
1-Abdominal Cramps
Press ↓ / ↑ Start +/- 100

Here we have chosen 100 steps at a time. Press up or down arrow to select program.



=== PROGRAMS ===
1-Abdominal Cramps
Press ↓ / ↑ Start ok

To start selected program press the start button until ok appears,
then press the up or down arrow to start the program.

Congratulations on your purchase of the Bioresonance Rife Frequency therapy Machine.

Now more powerful than ever. You can select the output voltage, (Between 5V to 26V-DC).

Including 2 copper rods. The Copper rods Comes with 3.5mm mono contact cable. And the copper rods is short-circuit protected.

**Do not open The Therapy machine as it contains sensitive electronics.
If you open it, the warranty no longer applies.**

For more information and spare parts visit www.bioresonance.nu

What is a Bioresonance Rife Frequency therapy machine?

This therapy is based on the principles of using specific frequencies to target and eliminate pathogens, such as bacteria, viruses, and parasites, within the body. The concept is that every organism has its own unique frequency, and by using the correct frequencies, these harmful organisms can be destroyed without harming healthy cells.

When combined in a single machine, Bioresonance Rife Frequency Therapy aims to provide a comprehensive approach to wellness by addressing both energetic imbalances within the body and targeting harmful pathogens. The goal is to support the body's natural healing mechanisms and promote overall health and vitality. Unlock Your Body's Healing Potential with our Bioresonance Rife Machine! Harnessing cutting-edge technology to target and eliminate imbalances at the cellular level. Experience holistic wellness and vitality like never before. Embrace the future of health today!

435 Pre-programmed "disease" programs with RIFE Frequences
(So it contains many 1000s of different frequencies).
Each program contains one to ten different frequencies.

All programs are according to Royal Raymond Rife's frequency list.

If a frequency program makes you feel uncomfortable, please stop running it until you feel better. Please note that these frequencies an the Machine are only for experiments.

Included in the package should be:

- 1: The bioresonator has a built-in rechargeable Litium battery of 10000mA 37W.
- 2: Biroresonance Machine. You can select the output voltage, Between 5V to 26V-DC.
- 3: 2 Copper rods with 3.5mm mono contact cable. The copper rods is short-circuit protected.
- 4: Charging cable USB-C.
- 5: Instruction manual.
- 6: Frequency list.